PERFORMANCE DATA SHEET

2610

\\nce||

Monoclonal anti-human Kappa Light chain*

mAb name/Clone: L1C1 *Isotype*: Mouse IgG1κ

Immunogen: Human B lymphoma cells

CATALOG#: 136-820 (Preservative-free)

QUANTITY: 100 μg CONCENTRATION: 1.0 mg/ml

INFORMATION: Human immunoglobulins are glycoproteins composed of two disulfide-bonded heavy (H) chain subunits, each of which is linked by interchain disulfide bonds to a light (L) chain forming a tetramolecular complex. There are five classes of immunoglobulins, designated IgG, IgA, IgM, IgD and IgE, which are defined by differences in the constant region of H chains. L chains are divided into kappa or lambda classifications based on structural antigenic differences. All classes of immunoglobulins have been found on the cell surface of B lymphocytes where they function as antigen receptors to elicit antigen-dependent proliferation and secretion of antigen specific soluble circulating antibodies. Antibody L1C1 recognizes cell surface kappa on B cells.

References: 1. Basic and Clinical Immunology, Seventh edition (D. P. Sites & A. I. Terr, eds.) Appleton & Lange., Norwalk, CT (1991).

STORAGE CONDITIONS: *Store at 2 - 5^oC*. Freeze/Thawing is not recommended. **Open under aseptic conditions.**

PRODUCT STABILITY: Product should retain activity for at least 12 months after shipping date when stored as recommended. Ship Date:_____


BUFFER: 50 mM Sodium Phosphate pH 7.5, 100 mM Potassium Chloride, 150mM NaCl.

PRODUCTION: Antibody was Protein A purified from (low FBS containing) tissue culture supernatant. Purity

was >95% Immunoglobulin by SDS-PAGE with less than 1% Bovine Immunoglobulin. Product was $0.2~\mu m$ filtered and vialed under aseptic conditions.

PERFORMANCE: Five x 10^5 ficoll prepared **human peripheral white blood cells** were pre incubated ~5 minutes with 20 µl of 250 µg/ml human Ig after which they were incubated 45 minutes on ice with 80 µl of anti-Kappa antibody at **5 µg/ml**. Cells were washed twice and incubated with 2^0 reagent Goat anti-Mouse IgG/FITC (Catalog #232-011), after which they were washed three times, fixed and analyzed by FACS. A net **28%** sub population of the cells stained positive with a mean shift of **1.20** \log_{10} fluorescent units when compared to a Mouse IgG1 negative control (Catalog #278-010).

Binding of anti-Kappa mAb +GAM/FITC to human PBMC

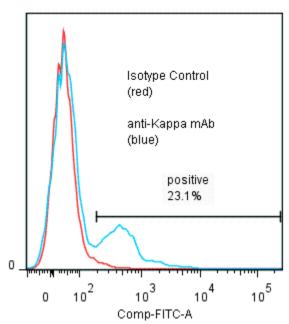

Ancell Corporation P.O. Box 87 Bayport, MN 55003-0087 USA Phone: Toll free 800-374-9523 or 651-439-0835 Fax: 651-439-1940

^{*} Research Use Only. Not for use in Diagnostic procedures.